

Алматы (7273)495-231	Иваново (4932)77-34-06	Магнитогорск (3519)55-03-13	Пермь (342)205-81-47	Тверь (4822)63-31-35
Ангарск (3955)60-70-56	Ижевск (3412)26-03-58	Москва (495)268-04-70	Ростов-на-Дону (863)308-18-15	Тольятти (8482)63-91-07
Архангельск (8182)63-90-72	Иркутск (395)279-98-46	Мурманск (8152)59-64-93	Рязань (4912)46-61-64	Томск (3822)98-41-53
Астрахань (8512)99-46-04	Казань (843)206-01-48	Набережные Челны (8552)20-53-41	Самара (846)206-03-16	Тула (4872)33-79-87
Барнаул (3852)73-04-60	Калининград (4012)72-03-81	Нижний Новгород (831)429-08-12	Саранск (8342)22-96-24	Тюмень (3452)66-21-18
Белгород (4722)40-23-64	Калуга (4842)92-23-67	Новокузнецк (3843)20-46-81	Санкт-Петербург (812)309-46-40	Ульяновск (8422)24-23-59
Благовещенск (4162)22-76-07	Кемерово (3842)65-04-62	Ноябрьск (3496)41-32-12	Саратов (845)249-38-78	Улан-Удэ (3012)59-97-51
Брянск (4832)59-03-52	Киров (8332)68-02-04	Новосибирск (383)227-86-73	Севастополь (8692)22-31-93	Уфа (347)229-48-12
Владивосток (423)249-28-31	Коломна (4966)23-41-49	Омск (3812)21-46-40	Симферополь (3652)67-13-56	Хабаровск (4212)92-98-04
Владикавказ (8672)28-90-48	Кострома (4942)77-07-48	Орел (4862)44-53-42	Смоленск (4812)29-41-54	Чебоксары (8352)28-53-07
Владимир (4922)49-43-18	Краснодар (861)203-40-90	Оренбург (3532)37-68-04	Сочи (862)225-72-31	Челябинск (351)202-03-61
Волгоград (844)278-03-48	Красноярск (391)204-63-61	Пенза (8412)22-31-16	Ставрополь (8652)20-65-13	Череповец (8202)49-02-64
Вологда (8172)26-41-59	Курск (4712)77-13-04	Петрозаводск (8142)55-98-37	Сургут (3462)77-98-35	Чита (3022)38-34-83
Воронеж (473)204-51-73	Курган (3522)50-90-47	Псков (8112)59-10-37	Сыктывкар (8212)25-95-17	Якутск (4112)23-90-97
Екатеринбург (343)384-55-89	Липецк (4742)52-20-81		Тамбов (4752)50-40-97	Ярославль (4852)69-52-93

Россия +7(495)268-04-70

Казахстан +7(7172)727-132

Киргизия +996(312)96-26-47

www.zoll.nt-rt.ru | | zof@nt-rt.ru

Технические характеристики на травматологические комплекты

Mobilize Trauma Kit

КОМПАНИИ **ZOLL**

Mobilize™ Trauma Kit

Comprehensive or Mobile kits



Box contents:

- Mobilize Trauma Kit
- Quick Start Guide (this document)
- Charging equipment (USB cable and wall adapter)
- Security seals and inspection card

Requirements:

Access to WiFi


Kit contents

A1	Tourniquet
B1	QuikClot®
B2	Pressure Dressing
C1	Chest Seal

D1	Face Shield
E1	Allergy Medicine
E2	Aspirin
E3	Glucose




1	Splint
2	Trauma Gauze
3	Gauze
4	Rolled Gauze
5	Biohazard Bag
6	Hand Wipes
7	Emergency Blanket
8	Ice Pack
9	Elastic Wrap
10	Burn Dressing
11	Triangular Bandage
PC	Portable Charger*

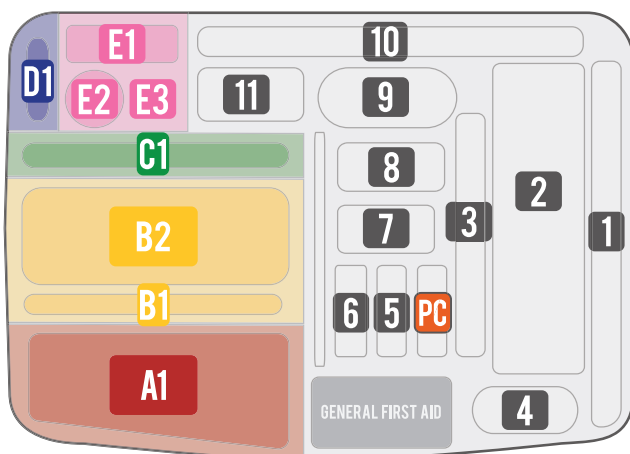*COMPREHENSIVE KIT ONLY



ON LID	
Gloves	Trauma Shears
Tape	Incident Pad & Pen

Setting up the included tablet

- 1) Turn the tablet on by pressing and holding the side button on the right until the tablet turns on (approximately 5 seconds).
- 2) After the Mobilize Rescue application launches, press "ACTIVATE APP."
- 3) Press the "SETTINGS"  button in the upper right-hand corner.
- 4) Press "Change WiFi" and toggle from Off to On.
- 5) Connect to your local WiFi network by selecting the network name and entering the password.
- 6) Press the back button  twice, then press the "ACTIVATE APP"  button to reach the Login/Register Page.
- 7) Press Register to create an account; or enter your existing account information, press "LOGIN" and skip to step 9.
- 8) Enter your new account information: First name, Last name, email, and password. The password must be at least 8 characters. Check the box "I have read and agreed with the Terms and conditions." Press "SUBMIT."
- 9) On the next page under "Activate app with," choose either:
 - a. License Code: provided by your local authorized ZOLL representative
 - b. Serial Number: Located near the hinge on the Comprehensive carrying case, and on the back of the tablet for Mobile units.
- 10) Press Submit and you will see "Take a Deep Breath, Mobilize Will Help You."
You are now ready to use Mobilize.

- 11) Use the charging equipment to charge the tablet's battery to 100%. To ensure that the kit is ready if needed at a moment's notice, charge the battery once per month.
 - The Comprehensive kit can be charged with the case closed. The tablet goes to sleep when lid is shut and wakes up when it is opened.
 - To put the Mobile kit to sleep, press the side button on the tablet. Double tap the screen to wake it up.

Please contact your Authorized ZOLL Distributor to purchase replenishments for utilized or expired supplies.

Public Access or Compact kits



Box contents:


Mobilize Trauma Kit
Quick Start Guide (this document)

Requirements:

Mobile phone or tablet for app installation
Access to WiFi or cellular connection

Kit contents

A1	Tourniquet	B2	Pressure Dressing
B1	QuikClot®	D1	Face Shield
C1	Chest Seal	7	Emergency Blanket

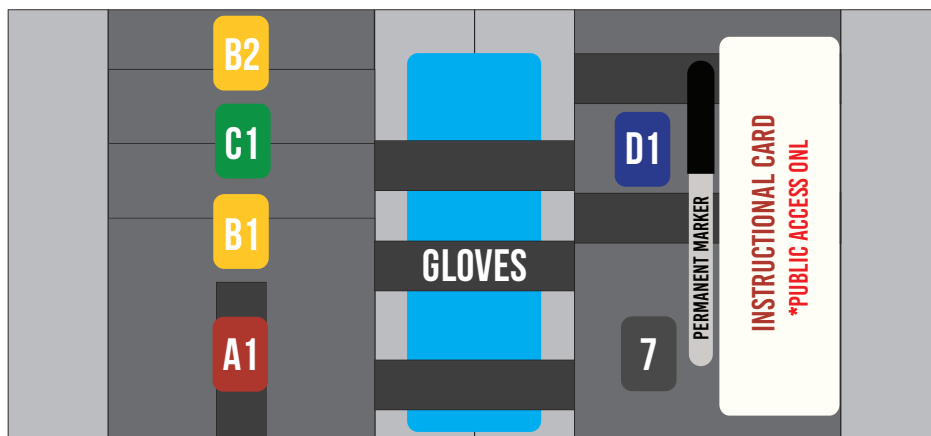

Setting up the mobile application

- 1) Go to the Apple App Store or Google Play Store and search for the Mobilize Rescue application.
- 2) Download the Mobilize Rescue application to your device.
- 3) When the download is complete, open the application.
(Apple users: if prompted to allow notifications, press "Allow.")
- 4) After Mobilize launches, press the "ACTIVATE APP" → button.
- 5) Press "REGISTER" to create an account; or enter your existing account information, press "LOGIN" and skip to step 7.
- 6) Enter your new account information: First name, Last name, email, and password. The password must be at least 8 characters. Check the box "I have read and agreed with the Terms and conditions." Press "SUBMIT."
- 7) On the next page under "Activate app with," choose either:
 - a. License Code: provided by your local authorized ZOLL representative
 - b. Serial Number: For the Public Access kit, the serial number is located on the front of the carrying case. For the Compact kit, use the number in the box here:
- 8) Press Submit and you will see "Take a Deep Breath, Mobilize Will Help You."
You are now ready to use Mobilize.

If you do not have a License Code, or wish to get access to the application for additional users, please contact your Authorized ZOLL Distributor.

Thank you for taking safety to heart and purchasing your Mobilize Trauma Kit (Mobile, Comprehensive, Public Access, or Compact). The Mobilize system is designed to provide lay rescuers with equipment and guidance to help save lives. Each kit includes medical and trauma supplies, and is paired with an interactive application that provides real-time instructions on how to use the supplies to provide care for injured or sick victims. The step-by-step directions provided through the Mobilize Rescue application enables bystanders to assess, monitor and manage medical emergencies in the critical minutes before EMS arrives.



Алматы (7273)495-231
Ангарск (3955)60-70-56
Архангельск (8182)63-90-72
Астрахань (8512)99-46-04
Барнаул (3852)73-04-60
Белгород (4722)40-23-64
Благовещенск (4162)22-76-07
Брянск (4832)59-03-52
Владивосток (423)249-28-31
Владикавказ (8672)28-90-48
Владимир (4922)49-43-18
Волгоград (844)278-03-48
Вологда (8172)26-41-59
Воронеж (473)204-51-73
Екатеринбург (343)384-55-89

Иваново (4932)77-34-06
Ижевск (3412)26-03-58
Иркутск (395)279-98-46
Казань (843)206-01-48
Калининград (4012)72-03-81
Калуга (4842)92-23-67
Кемерово (3842)65-04-62
Киров (8332)68-02-04
Коломна (4966)23-41-49
Кострома (4942)77-07-48
Краснодар (861)203-40-90
Красноярск (391)204-63-61
Курск (4712)77-13-04
Курган (3522)50-90-47
Липецк (4742)52-20-81

Магнитогорск (3519)55-03-13
Москва (495)268-04-70
Мурманск (8152)59-64-93
Набережные Челны (8552)20-53-41
Нижний Новгород (831)429-08-12
Новокузнецк (3843)20-46-81
Ноябрьск (3496)41-32-12
Новосибирск (383)227-86-73
Омск (3812)21-46-40
Орел (4862)44-53-42
Оренбург (3532)37-68-04
Пенза (8412)22-31-16
Петрозаводск (8142)55-98-37
Псков (8112)59-10-37

Пермь (342)205-81-47
Ростов-на-Дону (863)308-18-15
Рязань (4912)46-61-64
Самара (846)206-03-16
Саранск (8342)22-96-24
Санкт-Петербург (812)309-46-40
Саратов (845)249-38-78
Севастополь (8692)22-31-93
Симферополь (3652)67-13-56
Смоленск (4812)29-41-54
Сочи (862)225-72-31
Ставрополь (8652)20-65-13
Сургут (3462)77-98-35
Сыктывкар (8212)25-95-17
Тамбов (4752)50-40-97

Тверь (4822)63-31-35
Тольятти (8482)63-91-07
Томск (3822)98-41-53
Тула (4872)33-79-87
Тюмень (3452)66-21-18
Ульяновск (8422)24-23-59
Улан-Удэ (3012)59-97-51
Уфа (347)229-48-12
Хабаровск (4212)92-98-04
Чебоксары (8352)28-53-07
Челябинск (351)202-03-61
Череповец (8202)49-02-64
Чита (3022)38-34-83
Якутск (4112)23-90-97
Ярославль (4852)69-52-93

Россия +7(495)268-04-70

Казахстан +7(7172)727-132

Киргизия +996(312)96-26-47

www.zoll.nt-rt.ru | | zof@nt-rt.ru